


"When someone asks for an assistant, I say, 'Pull out a Thompson!'"

- Robert Mikkelsen, MD

Anterior Lumbar Articulating Arm System

Thompson's table-mounted Anterior Lumbar Articulating Arm Retractor System contains streamlined components for fast, easy set up. This specialized system utilizes easily adjustable articulating arms to attach retractor blades to a stable, table mounted frame. Articulating arms provide maximum positioning versatility, with the added benefit of multiplanar retraction.

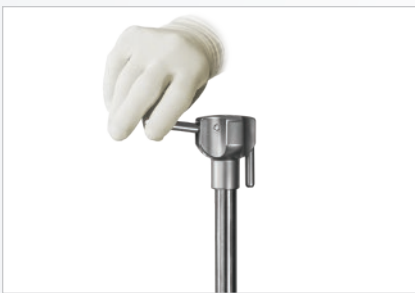


Sterile Mount

Top control allows attachment to table, as well as adjustments, to be achieved within the sterile field.



Fast, Easily Adjustable Arms


Quick, independent adjustments with maximum positioning versatility, including multiplanar retraction.

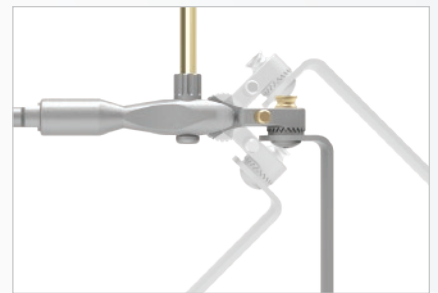

S-Lock® Quick Angle Technology

Allows quick, precise blade placement while remaining low profile and out of the way.

Anterior Lumbar Articulating Arm System ORDERING INFORMATION

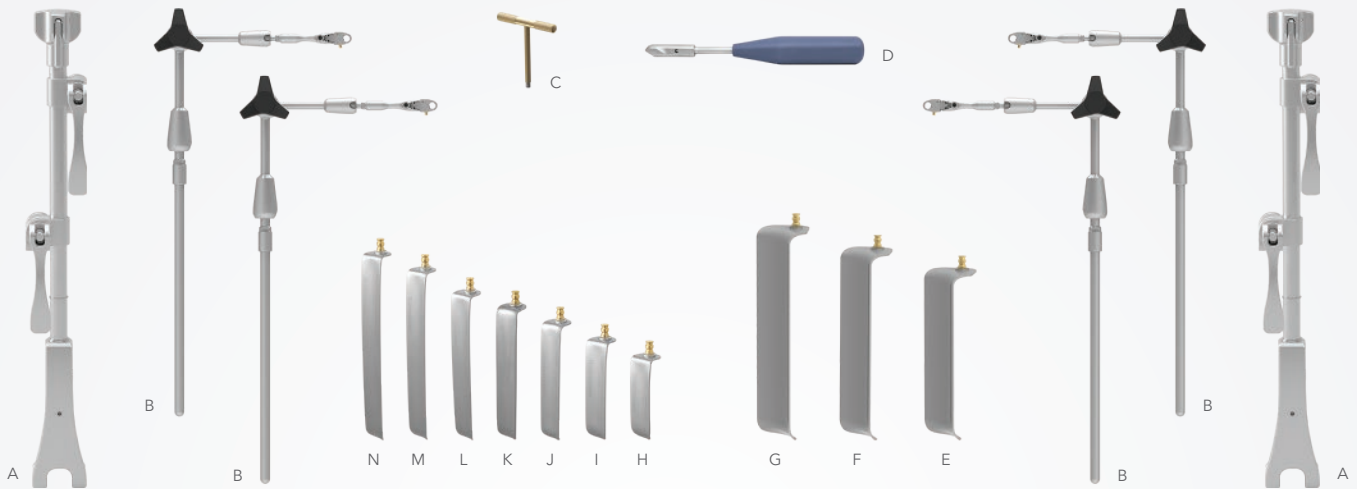

ANTERIOR LUMBAR ARTICULATING ARM SYSTEM #SL82048

REF	PART #	ITEM DESCRIPTION	QTY
Tray 1			
A	43903ACL	Elite III Rail Clamp with 2 Cam Joints 16"	1
B	SL42191G	Quick Angle Articulating Arm 15" x 15 1/2"	2
C	60020	T-Handle	1
D	42228	Manipulator Handle	1
	50000G	Instrument Case (23" x 11" x 3 1/2")	1
Tray 2			
A	43903ACL	Elite III Rail Clamp with 2 Cam Joints 16"	1
B	SL42191G	Quick Angle Articulating Arm 15" x 15 1/2"	2
C	60020	T-Handle	1
D	42228	Manipulator Handle	1
	50000G	Instrument Case (23" x 11" x 3 1/2")	1
E	SL46432ET	Radiolucent Richardson 50mm x 160mm (2" x 6 1/4")	1

REF	PART #	ITEM DESCRIPTION	QTY
F	SL46434ET	Radio. Richardson 50mm x 180mm (2" x 7")	1
G	SL46436ET	Radio. Richardson 50mm x 200mm (2" x 7 7/8")	1
Tray 3			
H	SL46558	Radiolucent Concave Short Neck Renal 25mm x 80mm (1" x 3 1/8")	4
I	SL46560	Radiolucent Concave Short Neck Renal 25mm x 100mm (1" x 4")	4
J	SL46562	Radiolucent Concave Short Neck Renal 25mm x 120mm (1" x 4 3/4")	4
K	SL46564	Radiolucent Concave Short Neck Renal 25mm x 140mm (1" x 5 1/2")	4
L	SL46566	Radiolucent Concave Short Neck Renal 25mm x 160mm (1" x 6 1/4")	4
M	SL46568	Radiolucent Concave Short Neck Renal 25mm x 180mm (1" x 7")	4
N	SL46570	Radiolucent Concave Short Neck Renal 25mm x 200mm (1" x 7 7/8")	4
	50000ALR	Instrument Case - AL Renal Vein Blades	1



Visit us online: thompsonsurgical.com
 For a Free Trial Call **1-800-227-7543** *

* Free trial valid for U.S. customers only. Customers outside the U.S.A., please call +1-231-922-0177 for availability.

© 2022 Thompson Surgical Instruments, Inc.

® The "T Circle" logomark is a Registered Trademark of Thompson Surgical Instruments, Inc.
 Patents: tpat.com

Rev A
 051622
 alaass

